

Endovascular aortic configurations

GORE® EXCLUDER® AAA Endoprosthesis

Trunk-Ipsilateral Leg Endoprosthesis

GORE® C3® Delivery System catalogue number	Intended aortic inner diameter (mm)	Aortic endoprosthesis diameter (mm)	Intended iliac inner diameter (mm)	Iliac endoprosthesis diameter (mm)	Overall device length (cm)	Recommended sheath size (Fr)
RLT231212	19-21	23	10-11	12	12	16
RLT231214	19-21	23	10-11	12	14	16
RLT231216	19-21	23	10-11	12	16	16
RLT231218	19-21	23	10-11	12	18	16
RLT231412	19-21	23	12-13.5	14.5	12	16
RLT231414	19-21	23	12-13.5	14.5	14	16
RLT231416	19-21	23	12-13.5	14.5	16	16
RLT231418	19-21	23	12-13.5	14.5	18	16
RLT261212	22-23	26	10-11	12	12	16
RLT261214	22-23	26	10-11	12	14	16
RLT261216	22-23	26	10-11	12	16	16
RLT261218	22-23	26	10-11	12	18	16
RLT261412	22-23	26	12-13.5	14.5	12	16
RLT261414	22-23	26	12-13.5	14.5	14	16
RLT261416	22-23	26	12-13.5	14.5	16	16
RLT261418	22-23	26	12-13.5	14.5	18	16
RLT281212	24-26	28.5	10-11	12	12	18
RLT281214	24-26	28.5	10-11	12	14	18
RLT281216	24-26	28.5	10-11	12	16	18
RLT281218	24-26	28.5	10-11	12	18	18
RLT281412	24-26	28.5	12-13.5	14.5	12	18
RLT281414	24-26	28.5	12-13.5	14.5	14	18
RLT281416	24-26	28.5	12-13.5	14.5	16	18
RLT281418	24-26	28.5	12-13.5	14.5	18	18
RLT311413	27-29	31	12-13.5	14.5	13	18
RLT311415	27-29	31	12-13.5	14.5	15	18
RLT311417	27-29	31	12-13.5	14.5	17	18
RLT351414	30-32	35	12-13.5	14.5	14	18
RLT351416	30-32	35	12-13.5	14.5	16	18
RLT351418	30-32	35	12-13.5	14.5	18	18

GORE® EXCLUDER® AAA Endoprosthesis

Contralateral Leg Endoprosthesis

Catalogue number	Intended iliac inner diameter (mm)	Iliac endoprosthesis diameter (mm)	Contralateral Leg length (cm)	Recommended sheath size (Fr)
PLC121000	10–11	12	10	12
PLC121200	10–11	12	12	12
PLC121400	10–11	12	14	12
PLC141000	12–13.5	14.5	10	12
PLC141200	12–13.5	14.5	12	12
PLC141400	12–13.5	14.5	14	12
PLC161000*	13.5–14.5	16	9.5	12
PLC161200*	13.5–14.5	16	11.5	12
PLC161400*	13.5–14.5	16	13.5	12
PLC181000*	14.5–16.5	18	9.5	12
PLC181200*	14.5–16.5	18	11.5	12
PLC181400*	14.5–16.5	18	13.5	12
PLC201000*	16.5–18.5	20	9.5	12
PLC201200*	16.5–18.5	20	11.5	12
PLC201400*	16.5–18.5	20	13.5	12
PLC231000*	18.5–21.5	23	10	14
PLC231200*	18.5–21.5	23	12	14
PLC231400*	18.5–21.5	23	14	14
PLC271000*	21.5–25.0	27	10	15
PLC271200*	21.5–25.0	27	12	15
PLC271400*	21.5–25.0	27	14	15

* PLEASE NOTE: All large diameter Contralateral Leg Endoprostheses (16, 18, 20, 23, 27 mm) can be used as Iliac Extenders.

GORE® EXCLUDER® AAA Endoprosthesis

Iliac Extender Endoprosthesis

Catalogue number	Intended iliac inner diameter (mm)	Endoprosthesis diameter (mm)	Endoprosthesis length (cm)	Recommended sheath size (Fr)
PLL161007	8–9	10	7	12
PLL161207	10–11	12	7	12
PLL161407	12–13.5	14.5	7	12

Aortic Extender Endoprosthesis

Catalogue number	Intended aortic inner diameter (mm)	Endoprosthesis diameter (mm)	Endoprosthesis length (cm)	Recommended sheath size (Fr)
PLA230300	19–21	23	3.3	16
PLA260300	22–23	26	3.3	16
PLA280300	24–26	28.5	3.3	16
PLA320400	27–29	32	4.5	17
PLA360400	30–32	36	4.5	18

GORE® EXCLUDER® AAA Endoprosthesis

Iliac Branch Component *

Catalogue number	Intended external iliac vessel diameter (mm)	Diameter iliac end (mm)	Iliac branch component length (cm)	Recommended sheath size (Fr)
CEB231010	6.5–9	10	10	16
CEB231210	10–11	12	10	16
CEB231410	12–13.5	14.5	10	16


Internal Iliac Component *

Catalogue number	Intended internal iliac vessel diameter (mm)	Diameter iliac end (mm)	Internal iliac component length (cm)	Recommended sheath size (Fr x cm)
HGB161007	6.5–9	10	7	12 x 45
HGB161207	10–11	12	7	12 x 45
HGB161407	12–13.5	14.5	7	12 x 45


* For use in the United States, add an 'A' at the end of the catalogue number (i.e., CEB231010A).

GORE® EXCLUDER® Iliac Branch Endoprosthesis


Deployment sequence

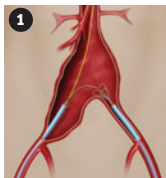


1 Establish sheaths and up-and-over throughwires via snare technique.

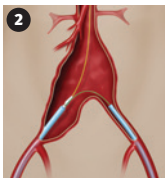


2 Advance Iliac Branch Component over both wires on GORE® EXCLUDER® Iliac Branch Endoprosthesis treatment side.

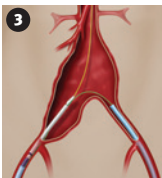


3 Withdraw sheath below GORE® EXCLUDER® Iliac Branch Endoprosthesis and align radiopaque marker ring above internal iliac ostium.

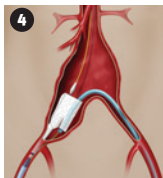


4 Deploy Iliac Branch Component (pull white outer deployment knob) and advance 12 Fr x 45 cm flexible sheath over bifurcation and into internal iliac gate.

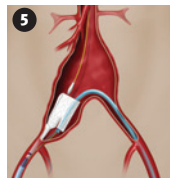


5 Cannulate internal iliac artery with third guidewire. Partially deployed Iliac Branch Component can be repositioned by: 1) rotating up to 90° in either direction, or 2) moving distally, to facilitate access to the internal iliac artery.

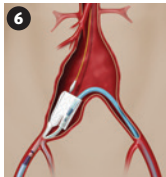


Advance Internal Iliac Component and align short Internal Iliac Component radiopaque marker with long Iliac Branch Component radiopaque marker.

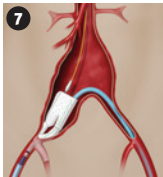


Deploy Internal Iliac Component and balloon overlap with Iliac Branch Component and seat within internal iliac artery. Remove up-and-over throughwire.

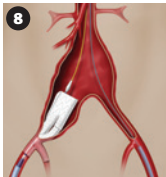

Deploy external iliac leg of Iliac Branch Component (pull gray inner deployment knob).

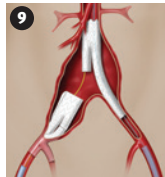

Exchange the 12 Fr flexible sheath for an 18 Fr sheath. Advance and deploy the Trunk-Ipsilateral Leg Endoprosthesis of the GORE® EXCLUDER® Device.

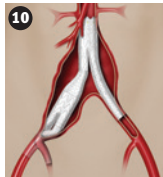

Advance and deploy 23 or 27 mm diameter Contralateral Leg as bridging component.

GORE® TAG® Conformable Thoracic Stent Graft with ACTIVE CONTROL System


Catalogue number	Intended aortic diameter (mm)	Proximal diameter (mm)	Distal diameter (mm)	Endoprosthesis length (cm)	Recommended GORE® Dryseal Flex Introducer Sheath size (Fr)	Partially uncovered stent length (mm)
TGM212110	16–19.5	21	21	10	18	3
TGM262110	19.5–24 / 16–19.5	26	21	10	20	4
TGM262610	19.5–24	26	26	10	20	4
TGM282810	22–26	28	28	10	20	4
TGM282815	22–26	28	28	15	20	4
TGMR312610	24–29 / 19.5–24	31	26	10	20	4
TGMR313110	24–29	31	31	10	20	4
TGMR313115	24–29	31	31	15	20	4
TGMR313120	24–29	31	31	20	20	4
TGM343410	27–32	34	34	10	22	5
TGM343415	27–32	34	34	15	22	5
TGM343420	27–32	34	34	20	22	5
TGMR373710	29–34	37	37	10	22	5
TGMR373715	29–34	37	37	15	22	5
TGMR373720	29–34	37	37	20	22	5
TGMR404010	31–37	40	40	10	22	6
TGMR404015	31–37	40	40	15	22	6
TGMR404020	31–37	40	40	20	22	6
TGM454510	34–42	45	45	10	24	6.5
TGM454515	34–42	45	45	15	24	6.5
TGM454520	34–42	45	45	20	24	6.5

Conformable GORE® TAG® Thoracic Endoprosthesis

Catalogue number	Intended aortic diameter (mm)	Proximal diameter (mm)	Distal diameter (mm)	Endoprosthesis length (cm)	Recommended sheath size (Fr)	Partially uncovered stent length (mm)
TGU212110	16–19.5	21	21	10	18	3
TGU262110	19.5–24 / 16–19.5	26	21	10	20	4
TGU262610	19.5–24	26	26	10	20	4
TGU282810	22–26	28	28	10	20	4
TGU282815	22–26	28	28	15	20	4
TGU312610	24–29 / 19.5–24	31	26	10	22	4
TGU313110	24–29	31	31	10	22	4
TGU313115	24–29	31	31	15	22	4
TGU343410	27–32	34	34	10	22	5
TGU343415	27–32	34	34	15	22	5
TGU343420	27–32	34	34	20	22	5
TGU373710	29–34	37	37	10	24	5
TGU373715	29–34	37	37	15	24	5
TGU373720	29–34	37	37	20	24	5
TGU404010	31–37	40	40	10	24	6
TGU404015	31–37	40	40	15	24	6
TGU404020	31–37	40	40	20	24	6
TGU454510	34–42	45	45	10	24	6.5
TGU454515	34–42	45	45	15	24	6.5
TGU454520	34–42	45	45	20	24	6.5

Intended range of aortic inner diameters (mm)

Device diameters

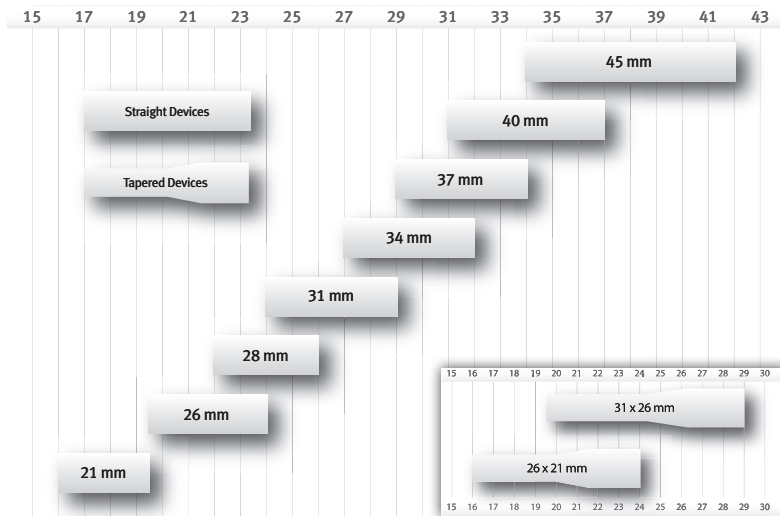

Aortic accessories

GORE® DrySeal Flex Introducer Sheath

Catalogue number	Sheath size (Fr)	Minimum sheath ID (mm)	Nominal sheath OD (mm)	Working length (cm)
DSF1233	12	4.0	4.7	33
DSF1245	12	4.0	4.7	45
DSF1433	14	4.7	5.3	33
DSF1533	15	5.0	5.6	33
DSF1633	16	5.3	6.1	33
DSF1833	18	6.0	6.7	33
DSF2033	20	6.7	7.5	33
DSF2065	20	6.7	7.5	65
DSF2233	22	7.3	8.2	33
DSF2265	22	7.3	8.2	65
DSF2433	24	8.0	8.8	33
DSF2465	24	8.0	8.8	65
DSF2633	26	8.7	9.5	33
DSF2665	26	8.7	9.5	65

Aortic accessories

GORE® Molding & Occlusion Balloon Catheter

Catalogue number	Size (Fr)	Balloon diameter (mm)	Catheter length (cm)
MOB37	10	10–37	90

GORE® Tri-Lobe Balloon Catheter

Catalogue number	Size (Fr)	Vessel diameter (mm)	Catheter length (cm)
BCM1634	18	16–34	104
BCL2645	18	26–42	104



Consult Instructions for Use

Products listed may not be available in all markets.

GORE®, ACTIVE CONTROL, C3®, EXCLUDER®, TAG®, and designs are trademarks of W. L. Gore & Associates. © 2017–2019 W. L. Gore & Associates, Inc.
AV1761-EN4 MAY 2019



W. L. GORE & ASSOCIATES, INC.
Flagstaff, AZ 86004

+65.67332882 (Asia Pacific)
00800.6334.4673 (Europe)
800.437.8181 (United States)
928.779.2771 (United States)

goremedical.com